

ZwickMaterials Testing

Product Information

Robotig Testing System 'roboTest L' (Linear) for Testing of Rubber O-Rings



Robotic testing system 'roboTest L' with testing machine 2.5 kN

Applications

The robotic testing system is used for the fully automatic performance of tensile tests on rubber rings (e.g. according to DIN 53504 with standard rings R1/R2, ISO 37 or ASTM D412).

System configuration

- Materials testing machine 2.5 kN up to 20 kN with roller grips according to DIN 53504 or ISO 37
- Removable specimen magazines for each 50 rings
- Robotic feeding system 'roboTest L' with pincer gripper
- Automatic thickness measuring unit for o-rings (measurement at one point)
- Industry Controller with test software testXpert® and automation software autoEdition2

Advantages of the Robotic Testing System 'roboTest L'


- A high reproducibility of the test results is obtained because operator influences are excluded (hand temperature, moist hands, eccentric or inclined insertion of specimens etc.).
- Qualified laboratory staff is relieved of routine jobs and is thus available for more complex activities.
- The machine can be used during idle times (break, night shift) and thus increases the rate of utilization and allows "quicker" results.
- The modular system makes an economical adaptation to specific customer requirements possible.
- The system reduces the testing costs per specimen and usually pays off within one to two years.
- Manual tests are still possible by simply moving the robotic feeding system aside.
- The usage of state-of-the-art web-technologies ensures a constant process control and remote diagnostics of the robotic testing system. Results as well as status messages can be sent directly per email or SMS.
- The automatic data logging system ensures secure documentation and enables statistical long-term monitoring (Statistical Process Control).
- The components of the robotic testing system are not subject to wear; they are maintenance-free and designed for three-shift operation.



ZwickMaterials Testing

Product Information

Robotig Testing System 'roboTest L' (Linear) for Testing of Rubber O-Rings

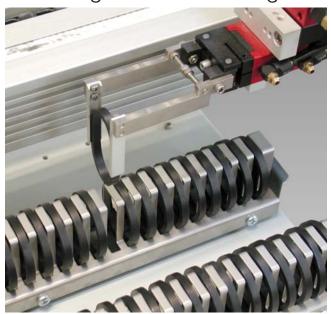

Pincer gripper takes O-ring out of magazine

Test Sequence

- The user fills the removable specimen magazine at the specimen preparation or directly on the test system with O-rings. Then the magazine will be placed on the magazine table.
- After the startup of the system, feed and test of the specimens are carried out automatically. With the optional expansion "motoric moving table" integrated, up to 7 specimen magazines (=294 specimens) are tested in succession.
- After testing of all magazined specimens the magazine can be refilled or changed by a prepared insert. A refilling of specimens in empty magazine places is possible at any time.
- By the individual loading of specimen into the magazines the order of testing can be controlled by the operator. Urgent specimens can be tested immediately by defining priorities.

Technical Data

Mechanics	
Mounting	coupled to the load frame
Dimensions (H x W x D)	1200 ¹⁾ x 700 ²⁾ x 1500 ³⁾ mm
Weight	approx. 200 kg (without
	specimens, depends on the
	equipment)

- 1) depends on the load frame
- 2) without option "motorized magazine table"
- 3) with linear axis

Connected values	
Electrical connection	230/115 V
Output	approx. 200 VA
Mains frequency	50/60 Hz
Compressed air	6 bar
Required compressed air	10 lpm

Control	
Automation	autoEdition2
Peripheral connection	RS 232

Specimens	
 Specimen type 	o-ring specimens
 Capacity 	50 (standard)
	294 (with option moving table)
 Material 	rubber
 Diameter (round) 	44.6 / 52.6 mm
 Thickness 	4 mm

Options

- Motorized magazine table for 294 specimens
- Data exchange with superior processor systems (e.g. LIMS) via upload/download of ASCII-files or ODBC
- Optical status indicator by threefold "traffic light" (running, refill specimens/finished, error)